

Temperature range*	18°C to 50°C
%RH *	50% to 90%
Temperature/%RH uniformity	1°C / 3%RH
% RH control at sensor	± 3%RH
Temperature control at sensor	± 1°C
Control resolution	± 0.1°C
Cooling/Heating type	Thermoelectric/Peltier driven-for heating/Cooling
Chamber Volume in Cu. Ft	43.2 cu.ft
External Dimension-LXWXH	51"×31"×81" (refer to drawing)
Internal Dimension-LXWXH	47" x 27" x 60" (refer to drawing)- height includes plenum which is 5"H.
Shelf dimension	23.5" x 24.5"
Cabinet specs	Stainless interior and exterior construction
Lights	Optional LED light optional
Electrical requirement	115V/60hz/1ph/6A, Unit plugs into standard wall outlet.
Shelf quantity, type	6 per chamber/ additional shelves available, wire shelf- PVC coated with 4 installation clips. Shelf capacity- 100 lbs distributed load.
Weight specs	Chamber weight- 340 lbs Shipping weight- 390 lbs

*Specifications are based on 20°C ambient and standard voltage. Specifications are subject to change without notice.

Chamber options:

- Chart recorder
- Additional shelving
- RO/DI water filtration system
- Water carboy

Chamber features:

Can den stability chamber is capable of maintaining stable temperature across wide operating temperature range of 10°C to 50°C and relative humidity range of 40% to 80% as limited by operational dew points and ambient condition. This chamber provides outstanding stability testing performance for manufacturers of pharmaceuticals, biologics, food, and personal care products. Can den stability chambers' performance meets ICH Q1A general case guidelines for pharmaceutical stability testing. Chamber is designed with plenum and a front to back airflow system for uniform sample exposure. Door locks are standard in doors. Chambers comes with optional 2" access port with insulation plug. Chamber uses thermoelectric assembly for heating and cooling. This type of cooling/heating system completely eliminates need for conventional refrigeration system and failures due to mechanical components such as compressors thereby drastically increasing the life expectancy of each chamber. Back wall plenum on chambers ensures uniform distribution of conditioned air. Chamber temperature set points are tightly controlled using PID controller. Rotronic humidity sensor or equivalent. Chamber is equipped with optional Future design touchscreen controller capable of showing live trends, data recorder.

Humidity system:




Humidity control is optional. Humidity range is achieved via additive ultrasonic humidification, a controller with with factory calibrated % RH sensor. Humidifier is installed typically in the rear of the chamber and can be easily serviced or replaced. RO/DI system is recommended but not required as long as the water quality specification is met. Tap water can be used if RO/DI water is not available. The humidifier comes with a 15' 1/2" polytube, shut off valve and quick disconnect fittings as well as a fitting to connect to water faucet.

Water quality specification:**Water quality specification:**

Water operating pressure range from 5 to 10psi


- Water Type: Deionized (DI) or Reverse Osmosis (RO) water.
- Target TDS Range: 15 ppm to 50 ppm.
- Conductivity Range: 30 µS/cm to 100 µS/cm.
- pH Range: Neutral (6.5 to 7.5).

Dehum: Dehumidification via customer facility supplied low dew point(-20c) compressed air that will hook up to our solenoid valve via 1/4" polytube at 20psi or optional Can den desiccant dehumidification blower system.

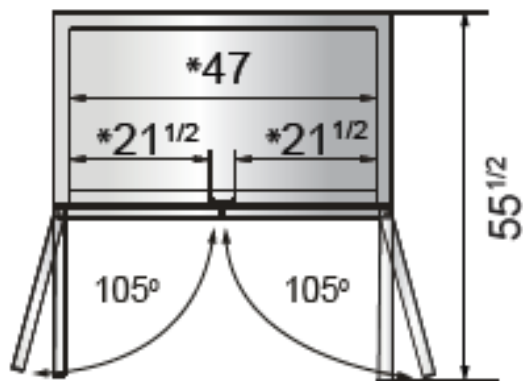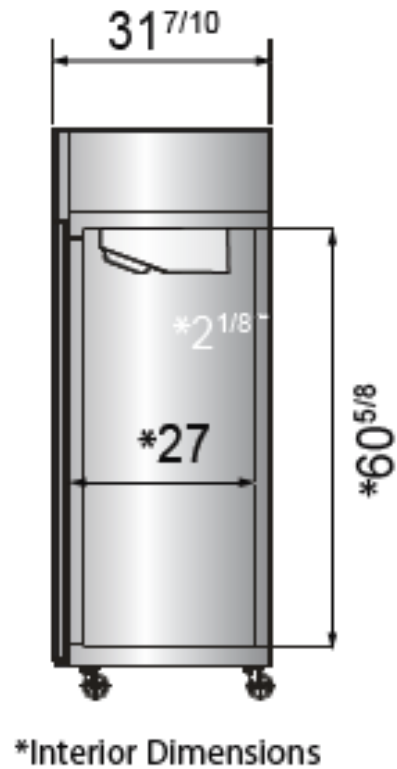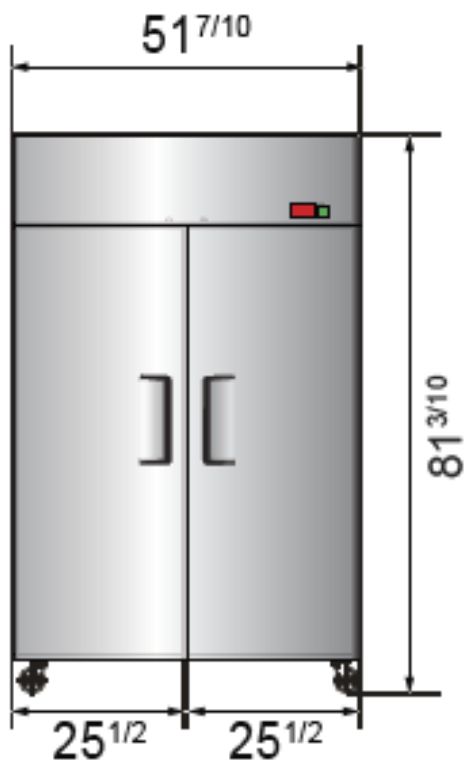

**SELF-CLOSING
AND STAY OPEN
DOOR FEATURE**

Facility RO or DI water valve

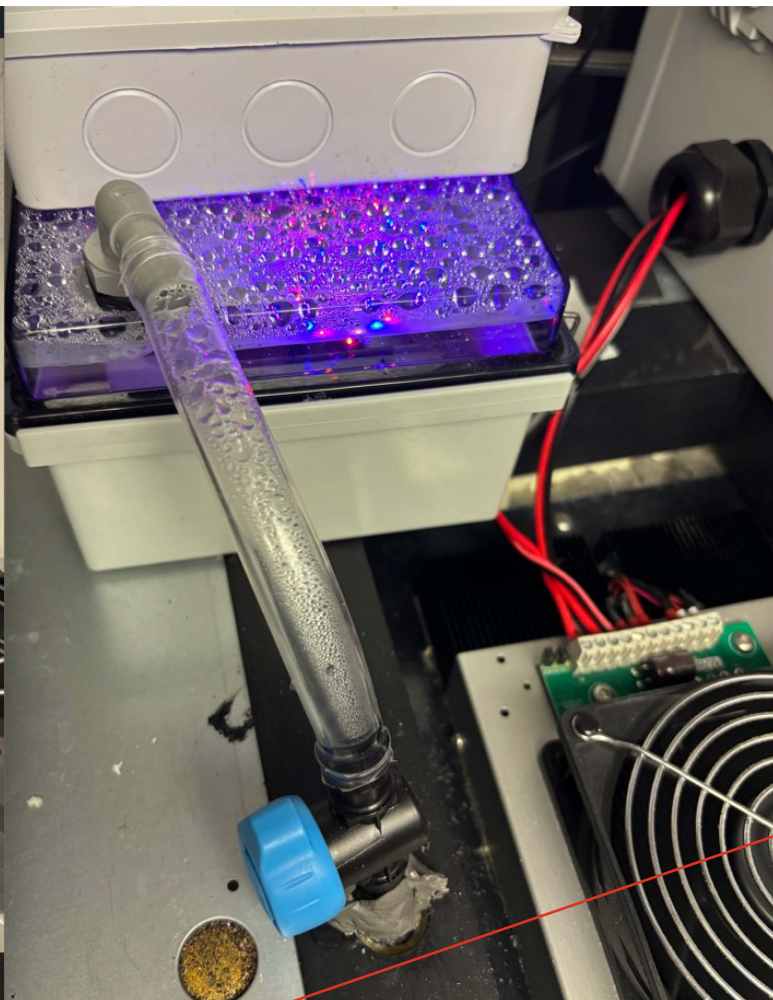

Facility RO or DI water line to humidifier--water pressure should be less than 10psi

Note: Connection line supplied by customer.

Connection from facility RO/DI water valve to humidifier done by customer.

Humidifier comes with 1/2" threaded adapter or 5/8" ID quick connect.

1/2 INCH FEMANLE NPT SIZE

